



imaging all the women



NEW

ICE
by RVF



The ultimate solution
to elevate your practice


delmont
imaging

Lead the way by setting new standards with **ICE** by RVF

Your unmet needs are now fully addressed!

Specifically designed for gynecology, ICE (innovative chips extractor) finally provides a fast and atraumatic solution for the extraction and recovery of resection chips. It redefines the standards of hysteroscopic resection by providing an effective response to a need that, until now, remained without a truly suitable solution.

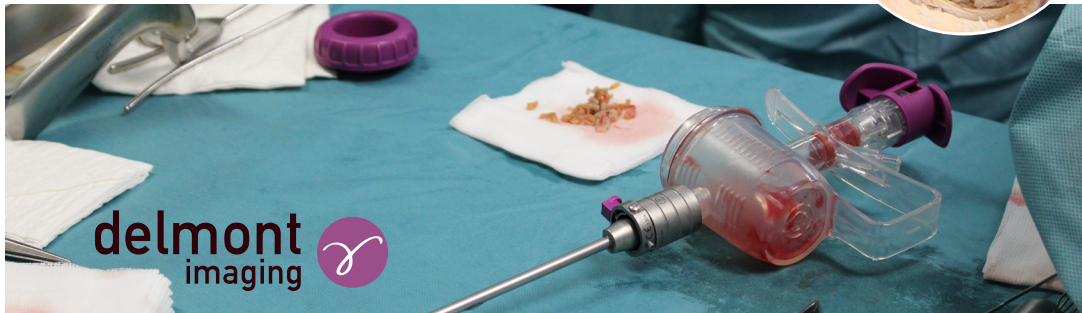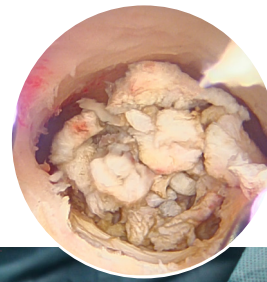

Innovation powered by patented technology

Created in close collaboration with experts in gynecology, ICE was designed by and for specialist surgeons. Guided by Professor Fernandez, whose initials sign this product, it meets requirements of ergonomics, ease of use and high efficiency.

"With just two presses of the hand, the cavity is completely emptied of chips but also blood, offering a considerable saving of time and an optimized operating day."


Prof. Hervé Fernandez.

The ideal approach for your patient care during resection

Retrieving chips at the end of surgery has always been a tedious and delicate task. Current techniques all have major issues: the use of the curette can lead to synechiae¹ while using the resectoscope itself may cause cervical lacerations. These methods will soon be a thing of the past!

Our ICE inner sheath, connected to the ICE activator, is fully compatible with the outer sheaths of ReseCare and Maxi-Care resectoscopes. This keeps the outer sheath in place inside the uterus throughout the complete procedure: the resection, the ICE installation and the chips extraction!

¹ Schenker JG. Intrauterine adhesions: an updated appraisal. *Fertil Steril.* mai 1982;37(5):593-610

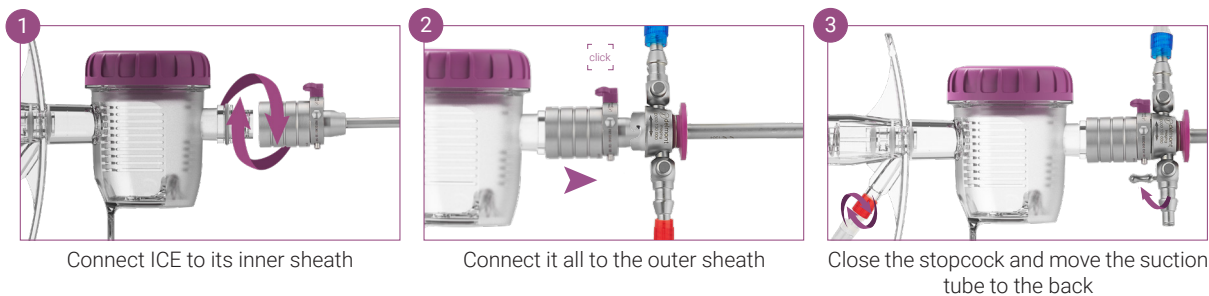

References :

ICE by RVF activator
Sterile. Single use. Packing by 5
Ref. D110 200 000

ICE by RVF inner sheath
For 26Fr. resectoscope
Ref. D300 130 107

ICE by RVF inner sheath
For 18.5Fr. resectoscope
Ref. D300 130 108



Zone Athélia IV - 297, Avenue du Mistral - 13600 La Ciotat - France
Tél : +33 (0) 9 51 51 30 30 - Fax : +33 (0) 9 57 51 31 00
www.delmont-imaging.com - contact@delmont-imaging.com



These medical devices are regulated health products that bear, under this regulation, the CE mark. They are manufactured by Delmont Imaging. They are intended for professional use only, read the instructions carefully before use. This document is not meant to replace the instruction for use.

Non contractual document. SAS Delmont Imaging, RCS MARSEILLE 817 710 411. All rights reserved. No information or part of this document may be reproduced or modified in any form without prior permission from Delmont Imaging.
REF D900 700 241 A-EN 2026-01.