

The complete **EasyCare** set – Ref. D400 100 000 includes the following instruments:

HD Hysteroscope, 2.9mm, 30°

Autoclavable. L:300mm

Ref. D300 100 000



Outer graduated operating sheath, 5.5mm

Rotatable suction stopcock. Quick lock. L: 216mm

Ref. D300 110 000



Inner graduated operating sheath, 4.9mm

7Fr. Operative canal, Rotatable irrigation stopcock. L: 230mm

Ref. D300 110 001



Graduated diagnostic sheath, 4mm

Rotatable irrigation stopcock. L: 270mm

Ref. D300 110 002



Graduated scissors, 7Fr.

Semi-rigid. L: 340mm. Single action jaw

Ref. D300 120 000



Graduated biopsy forceps, 7Fr.

Semi-rigid. L: 340mm. Double action jaws

Ref. D300 120 002



Graduated forceps (alligator), 7Fr.

Semi rigid. L: 340mm. Double action jaws

Ref. D300 120 001



Graduated IUD grasping forceps, 7Fr.

Semi rigid. L: 340mm. Double action jaws

Ref. D300 120 003



Light cable

Autoclavable. Ø:3.5mm. L: 2.3m. delivered with Storz connectors

Ref. D200 150 500



Optional accessories:

Bipolar electrode, flexible, button, 5Fr.

Ref. D300 140 010



Bipolar electrode, flexible, with tip, 5Fr.


Ref. D300 140 011



Sterilization basket

Autoclavable. 480x250x47mm. Supplied with a silicone mat and a lid


Ref. D200 160 500

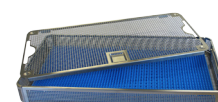

Zone Athélia V - 390, Avenue du Mistral - 13600 La Ciotat - France
Tel : +33 (0) 9 51 51 30 30 - Fax : +33 (0) 9 57 51 31 00
www.delmont-imaging.com - contact@delmont-imaging.com



imaging all the women

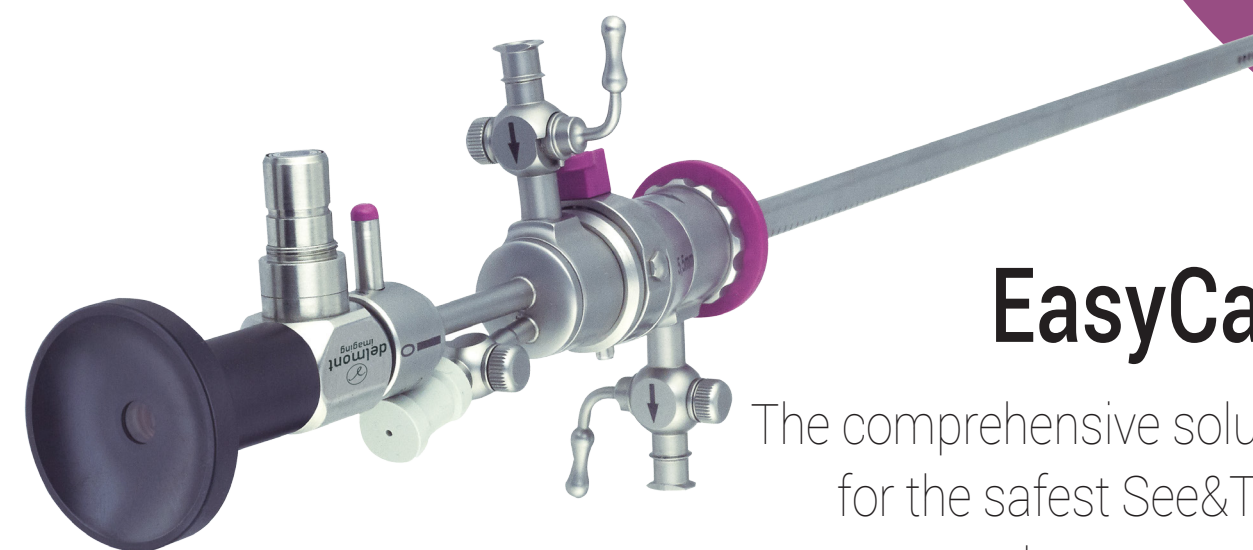

EasyCare

The comprehensive solution
for the safest See&Treat
hysteroscopies

According to Pr. H. Fernandez




Non contractual document - Ref. D900 700 010 C - 06/2018 - Copyright © 2016 DELMONT IMAGING. All rights reserved. No information or any part of this document may be reproduced or transmitted in any form without the prior permission of DELMONT IMAGING.



Pr. H. Fernandez,
Head of the Department of OB/GYN
of the Kremlin Bicêtre Hospital, France

« Many procedures can be performed in consulting room during a diagnostic hysteroscopy, avoiding your patients to come back for hospitalization under anesthesia. This is one of the many benefits of the See&Treat protocole, as well as the immediate care of the patient and the reduction of operating theatre programs, making room for heavier surgeries. To achieve See&Treat serenely, it is necessary to have a comprehensive and efficient equipment. **EasyCare** offers the right solution with smart sheaths but above all with a set of 7Fr. instrumentation much more efficient than the traditional 5Fr. one. Your gesture is precise with an immediate result and your patient satisfied. In addition, each instrument is graduated to ease your surgical marks. »

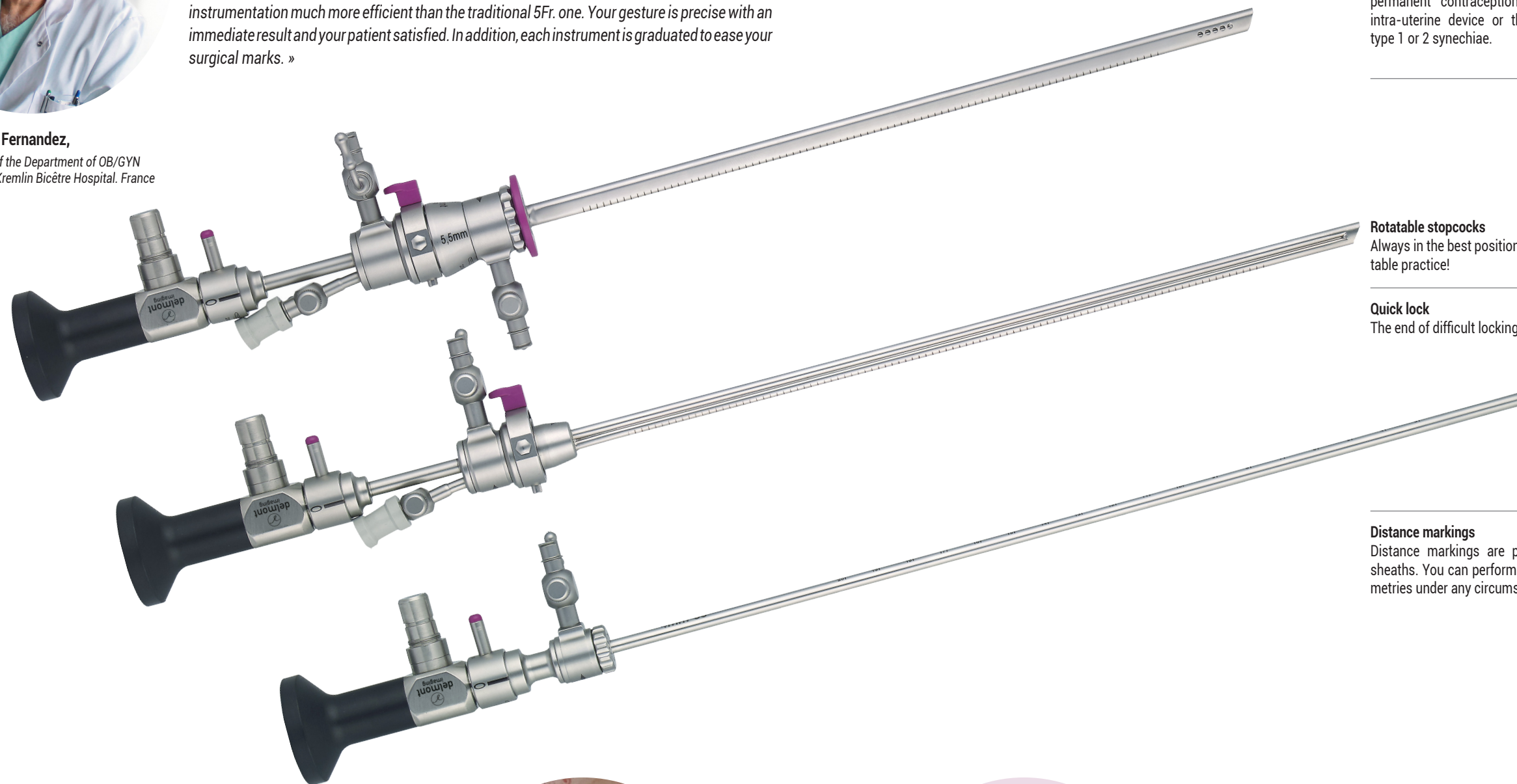


Inner atraumatic 4.9 mm sheath

With its blunt tip and even more restrained diameter, the inner sheath has been designed to be used without the outer sheath. It is the ideal tool for acts requiring no aspiration such as permanent contraception, removal of intra-uterine device or the section of type 1 or 2 synechiae.



Rotatable stopcocks

Always in the best position for a comfortable practice!

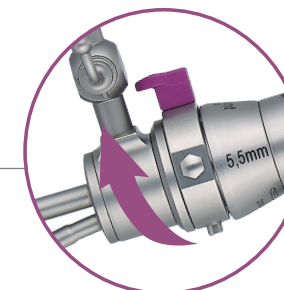


Quick lock

The end of difficult lockings.



Distance markings

Distance markings are present on all sheaths. You can perform your hysteroscopies under any circumstances.

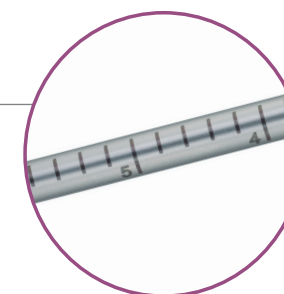

HD Hysteroscope : Unmatched performances and outstanding durability

Our latest generation of lenses, designed to take full advantage of HD, provides detailed and perfectly clear images even in the backgrounds. This hyper resolution together with a real homogeneous light distribution provides a huge depth of field. The whole image will be sharp! Our HD hysteroscope is of course autoclavable. Its three separate tubes design makes it resistant, durable but also easily repairable in case of accident.

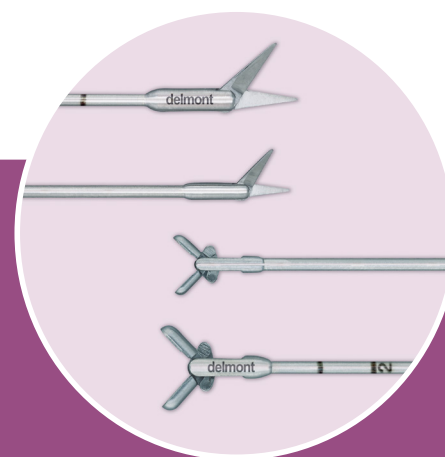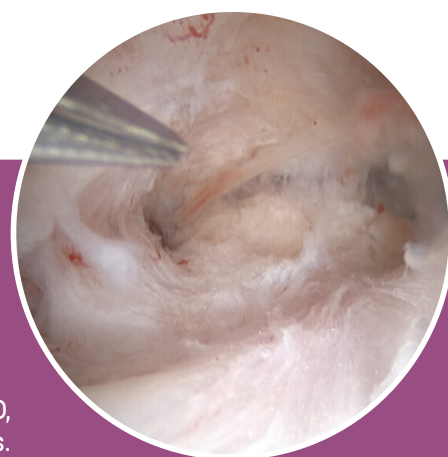

The 7Fr.¹ efficiency!

The 7Fr. scissors blades are 50% longer² than the 5Fr. ones, and the volume of material collected with a 7Fr. biopsy forceps is 3 times higher³ than the 5Fr. !

The 7Fr. instrumentation is naturally more rigid and more precise than the 5Fr. Here again, you save time and keep your patients more comfortable.

¹: 7mm = 3Fr.

²: 5.7mm vs 3.8mm

³: 8.23mm³ vs 3.2mm³