

The complete **MiniCare** set - Ref. D400 100 045 includes the following products:

**MiniCare operative hysteroscope. 6mm. 12°**

45° offset eyepiece. Autoclavable  
10Fr. working channel

Ref. D300 100 049



**18.5Fr. outer sheath**

Double flow. Quick-lock

Ref. D300 130 000



**Plastic handle for 3mm instruments**

With rotation module. Without ratchet.  
4 items are included in the set

Ref. D300 120 533



**3mm straight micro-scissors**

Single action jaws. Working length: 330mm

Ref. D300 120 509



**3mm curved scissors. Metzenbaum type**

Double action jaws. Working length: 330mm

Ref. D300 120 507



**3mm universal grasping forceps**

Toothed and wavy jaws  
Double action jaws. Working length: 330mm

Ref. D300 120 510



**3mm claw grasping and extracting forceps. With 2/3 teeth**

Single action jaws. Working length: 330mm

Ref. D300 120 515



**Light cable**

Autoclavable. Ø:3,5mm. L:2,3m. Delivered with Storz connectors

Ref. D200 150 500



**MiniCare sterilisation basket**

Includes two associable trays and a lid.  
With personalised silicone holders

Ref. D200 160 521



Optional accessories:

**Ø 3mm L-Hook electrode**

Working length: 350mm

Ref. D300 140 133



**Ø 3mm straight electrode**

Working length: 350mm

Ref. D300 140 134



**Ø 3mm V-Hook electrode**

Working length: 350mm

Ref. D300 140 135





*imaging all the women*



## MiniCare

The product that had to be invented  
for patients of childbearing age





## The perfect match between preservation and performance


To treat all uterine pathologies without compromise

MiniCare will seduce you for the treatment of your patients, especially those of childbearing age. Its small outer diameter reduces dilation steps and eases insertion, even with nulliparous patients.

The handling is easy and your surgery accurate.

The use of 3mm instruments without energy preserves the endometrium of your patients. Even in case of major pathologies such as large calcified fibroids or complete uterine septa, the performances of these rigid instruments allow you quick and efficient treatment.

**Offset optics**  
For a straight working channel  
allowing the use of rigid instruments  
while maintaining HD vision quality!

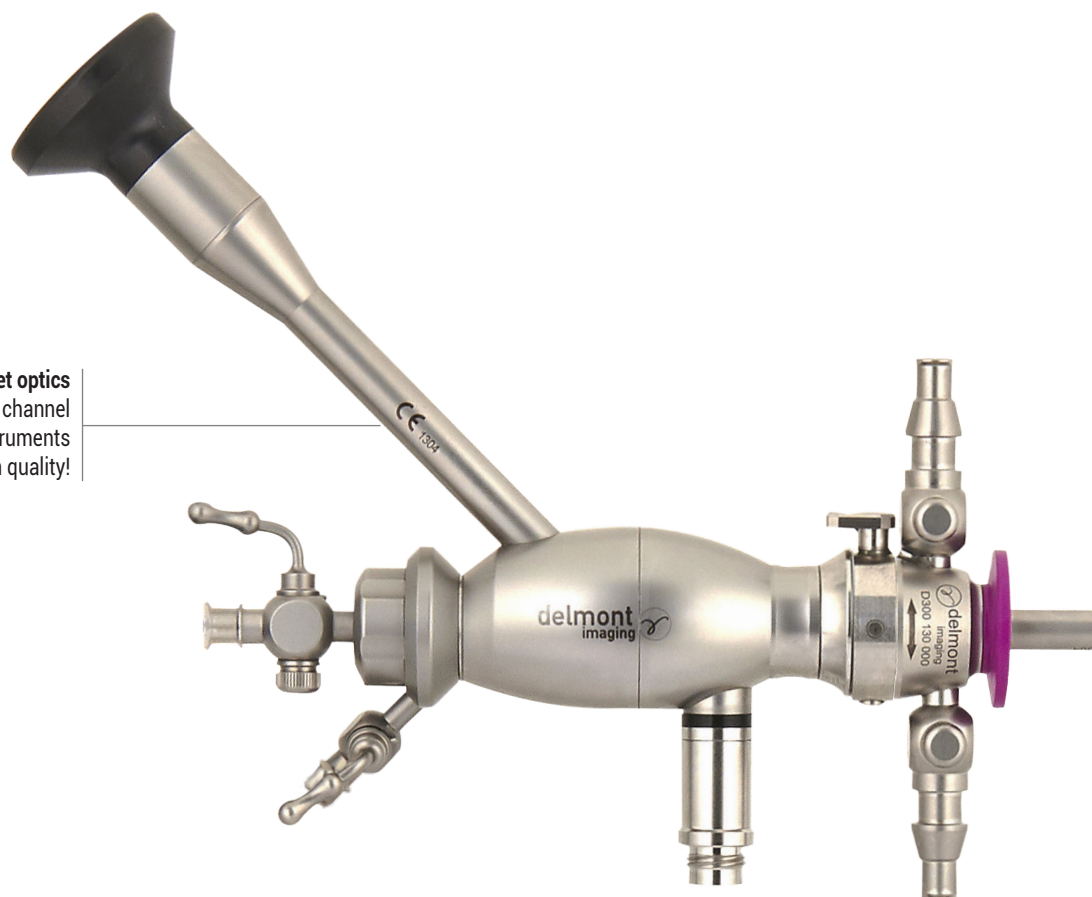

## Technical requirements pushed down to the smallest detail

The MiniCare operating hysteroscope has been designed to ensure the best performances and ergonomics.

The unique design of our hysteroscope gives you a natural grip, freeing up your workspace for better instruments handling.

The combination of HD lenses with generous and homogeneous lighting guarantees you a bright HD image.




**Three shapes of bipolar electrodes available**

Keep the safety of bipolar coagulation at your fingertips. These new generation electrodes can also be used with a cutting current.




**One single operating time, even for complex surgeries**

Find the power of your mini-laparoscopy instruments during your hysteroscopy. The efficiency and precision of rigid 3mm instruments reduce the duration of procedures.

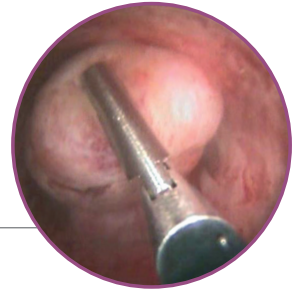

**Reduced diameter and safety**

The small external diameter is the assurance of less dilation of the cervix, minimizing the risk of false paths.

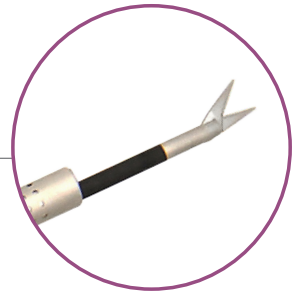

**MiniCare+<sup>1</sup>, the 2-in-1 instrumentation set!**

The outer sheath of MiniCare is identical to that of our bipolar resectoscope ReseCare. These two products being complementary, you can combine their use during the same surgery. The outer sheath remains positioned in the uterus during the instrument change, avoiding any cervical lacerations.

<sup>1</sup> : MiniCare+ = MiniCare and ReseCare