



imaging all the women



MiniCare

The product that had to be invented
for patients of childbearing age




The perfect match between preservation and performance

To treat all uterine pathologies without compromise

MiniCare will seduce you for the treatment of your patients, especially those of childbearing age. Its small outer diameter reduces dilation steps and eases insertion, even with nulliparous patients. The handling is easy and your surgery accurate.

The use of 3mm instruments without energy preserves the endometrium of your patients. Even in case of major pathologies such as large calcified fibroids or complete uterine septa, the performances of these rigid instruments allow you quick and efficient treatment.

Offset optics
For a straight working channel allowing the use of rigid instruments while maintaining HD vision quality!

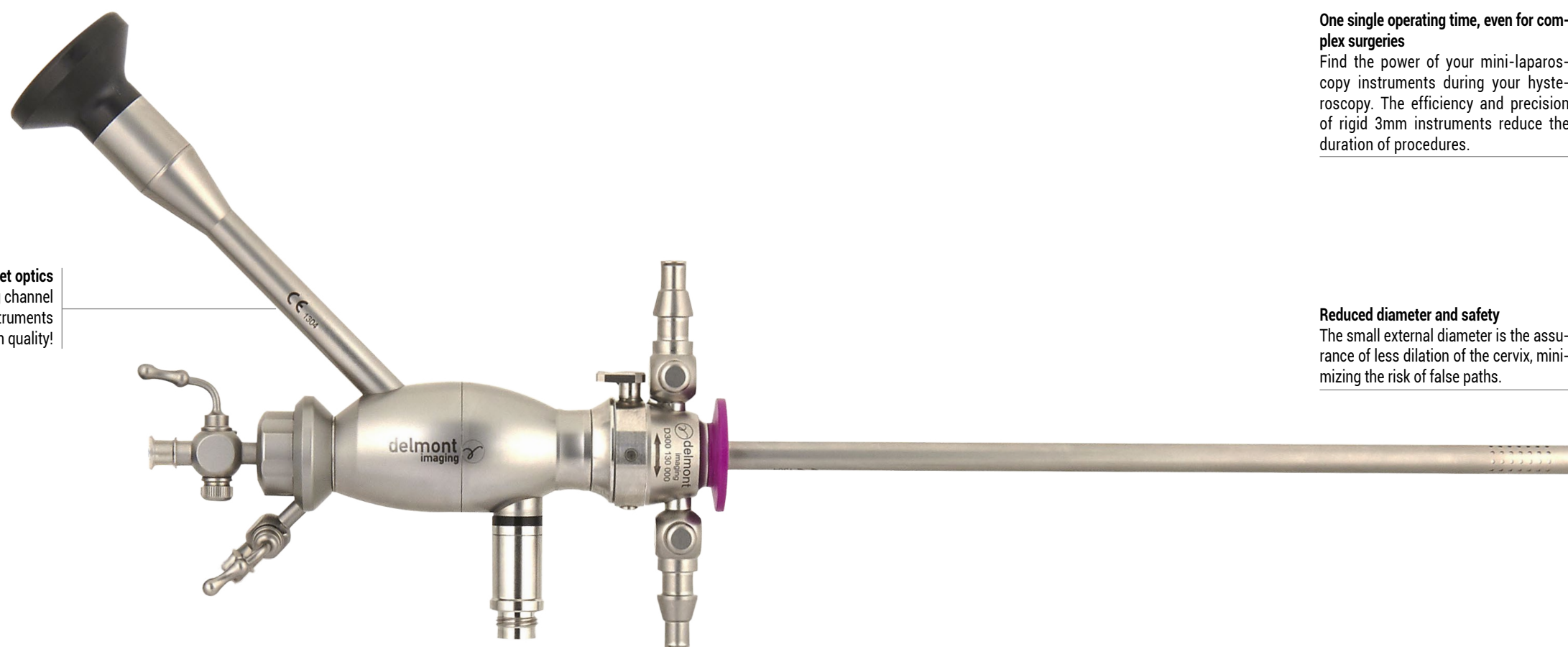


Three shapes of bipolar electrodes available

Keep the safety of bipolar coagulation at your fingertips. These new generation electrodes can also be used with a cutting current.




One single operating time, even for complex surgeries

Find the power of your mini-laparoscopy instruments during your hysteroscopy. The efficiency and precision of rigid 3mm instruments reduce the duration of procedures.

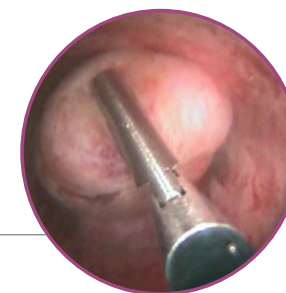

Reduced diameter and safety

The small external diameter is the assurance of less dilation of the cervix, minimizing the risk of false paths.

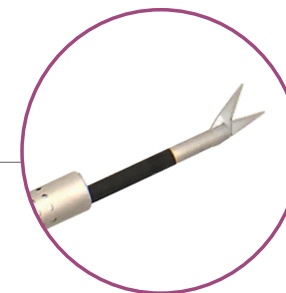

Technical requirements pushed down to the smallest detail

The MiniCare operating hysteroscope has been designed to ensure the best performances and ergonomics.

The unique design of our hysteroscope gives you a natural grip, freeing up your workspace for better instruments handling.

The combination of HD lenses with generous and homogeneous lighting guarantees you a bright HD image.



MiniCare+¹, the 2-in-1 instrumentation set!

The outer sheath of MiniCare is identical to that of our bipolar resectoscope ReseCare. These two products being complementary, you can combine their use during the same surgery. The outer sheath remains positioned in the uterus during the instrument change, avoiding any cervical lacerations.

¹ : MiniCare+ = MiniCare and ReseCare

The complete **MiniCare** set - Ref. D400 100 045 includes the following products:

MiniCare operative hysteroscope. 6mm. 12°

45° offset eyepiece. Autoclavable
10Fr. working channel

Ref. D300 100 049



18.5Fr. outer sheath

Double flow. Quick-lock

Ref. D300 130 000



Plastic handle for 3mm instruments

With rotation module. Without ratchet.
4 items are included in the set

Ref. D300 120 533



3mm straight micro-scissors

Single action jaws. Working length: 330mm

Ref. D300 120 509



3mm curved scissors. Metzenbaum type

Double action jaws. Working length: 330mm

Ref. D300 120 507



3mm universal grasping forceps

Toothed and wavy jaws
Double action jaws. Working length: 330mm

Ref. D300 120 510



3mm claw grasping and extracting forceps. With 2/3 teeth

Single action jaws. Working length: 330mm

Ref. D300 120 515



Light cable

Autoclavable. Ø:3,5mm. L:2,3m. Delivered with Storz connectors

Ref. D200 150 500



MiniCare sterilisation basket

Includes two associable trays and a lid.
With personalised silicone holders

Ref. D200 160 521



Optional accessories:

Ø 3mm L-Hook electrode

Working length: 350mm

Ref. D300 140 133



Ø 3mm straight electrode

Working length: 350mm

Ref. D300 140 134



Ø 3mm V-Hook electrode

Working length: 350mm

Ref. D300 140 135

