




*imaging all the women*

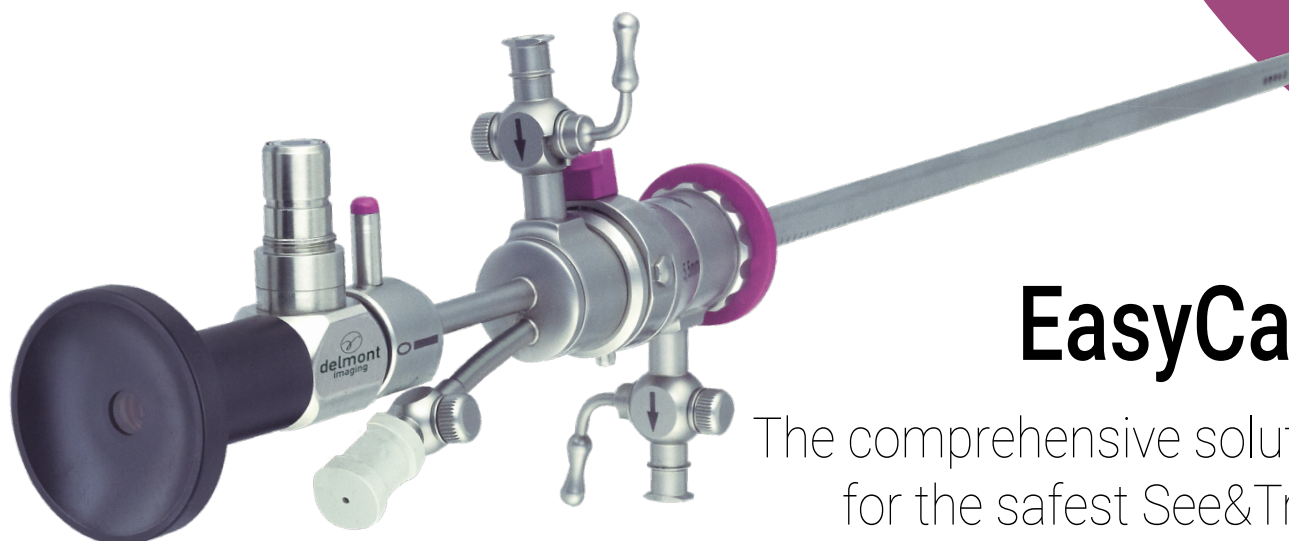

## EasyCare

The comprehensive solution  
for the safest See&Treat  
hysteroscopies


*According to Pr. H. Fernandez*



« Many procedures can be performed in consulting room during a diagnostic hysteroscopy, avoiding your patients to come back for hospitalization under anesthesia. This is one of the many benefits of the See&Treat protocole, as well as the immediate care of the patient and the reduction of operating theatre programs, making room for heavier surgeries. To achieve See&Treat serenely, it is necessary to have a comprehensive and efficient equipment. **EasyCare** offers the right solution with smart sheaths but above all with a set of 7Fr. instrumentation much more efficient than the traditional 5Fr. one. Your gesture is precise with an immediate result and your patient satisfied.

**Pr. Hervé Fernandez.** President of Delmont imaging scientific council. World-renowned hysteroscopy expert.

## EasyCare. One set, three possible tools.

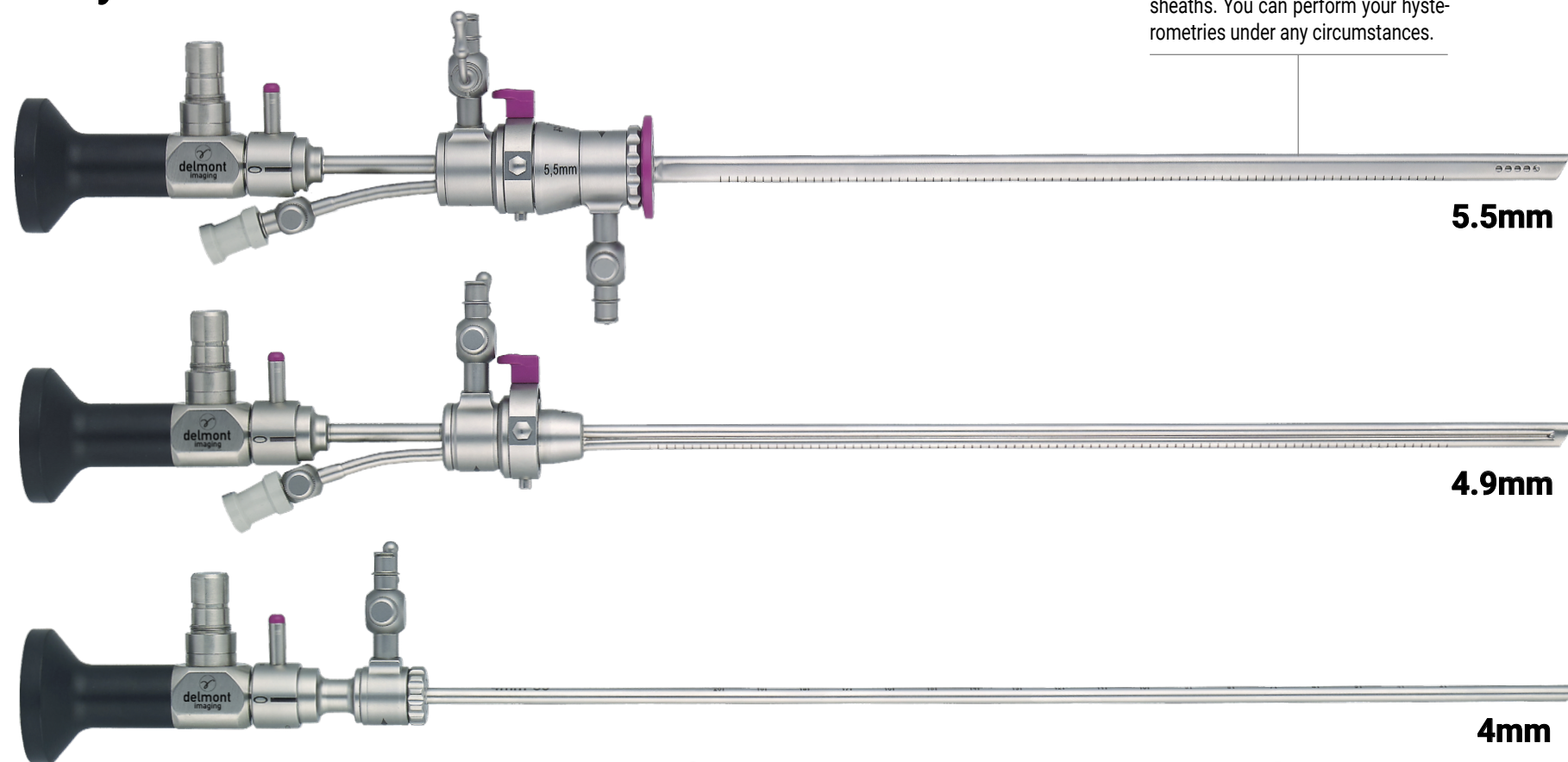


### Distance markings

Distance markings are present on all sheaths. You can perform your hysteroscopies under any circumstances.


### Inner atraumatic 4.9 mm sheath

With its blunt tip and even more restrained diameter, the inner sheath has been designed to be used without the outer sheath. It is the ideal tool for acts requiring no aspiration such as permanent contraception, removal of intra-uterine device or the section of type 1 or 2 synechiae.

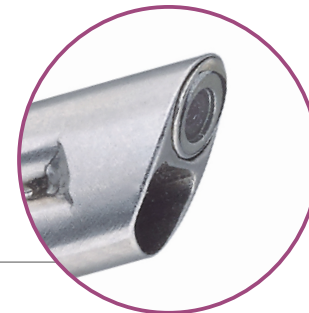

### Rotatable stopcocks

To never be hampered and always find the best position, direct each stopcock independently!

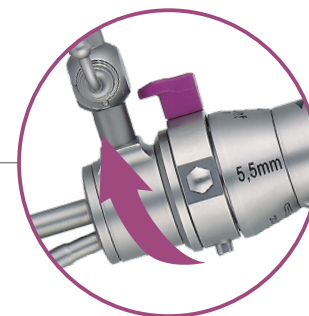


### Quick lock

The end of difficult lockings.

### Rotative semi-rigid instruments

Our instruments include a new rotation module facilitating use and greatly reducing the risk of torsion. Your comfort and the lifetime of your instruments are thus improved.


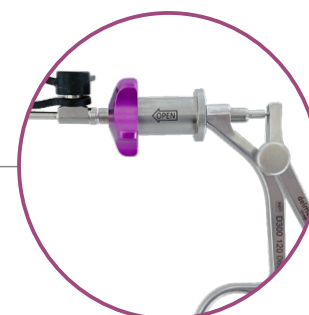


### Thought out down to the last detail

All use phases of the product are simplified: assembly and disassembly of the handle in one movement via the quick lock, as well as cleaning thanks to the easily accessed Luer-lock channel.

## HD Hysteroscope : Unmatched performances and outstanding durability

Our latest generation of lenses, designed to take full advantage of HD, provides detailed and perfectly clear images even in the backgrounds. This hyper resolution together with a real homogeneous light distribution provides a huge depth of field. The whole image will be sharp! Our HD hysteroscope is of course autoclavable. Its three separate tubes design makes it resistant, durable but also easily repairable in case of accident.

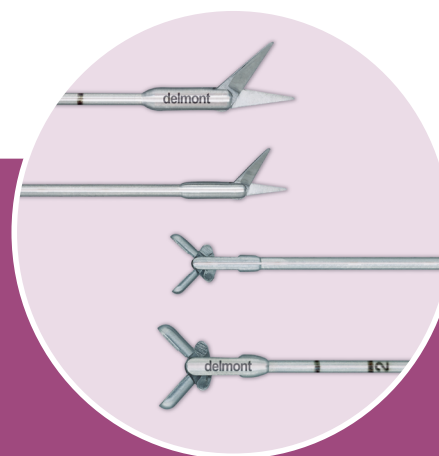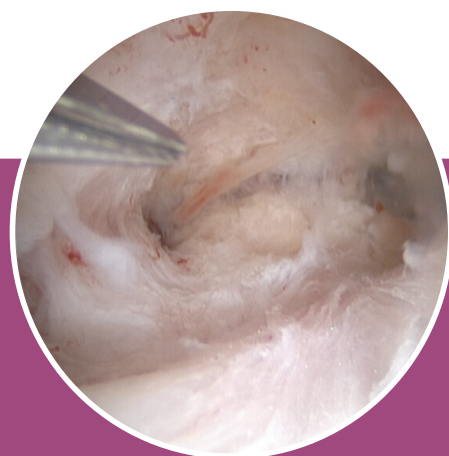

## The 7Fr.<sup>1</sup> efficiency!

The 7Fr. scissors blades are 50% longer<sup>2</sup> than the 5Fr. ones, and the volume of material collected with a 7Fr. biopsy forceps is 3 times higher<sup>3</sup> than the 5Fr. !

The 7Fr. instrumentation is naturally more rigid and more precise than the 5Fr. Here again, you save time and keep your patients more comfortable.

<sup>1</sup>: 1mm = 3Fr.

<sup>2</sup>: 5.7mm vs 3.8mm

<sup>3</sup>: 8.23mm<sup>3</sup> vs 3.2mm<sup>3</sup>

# The complete **EasyCare+** set – Ref. D400 100 017 includes the following instruments:

## HD Hysteroscope. 2.9mm. 30°

Autoclavable. L:300mm

Ref. D300 100 000



## Graduated diagnostic sheath. 4mm

Rotatable irrigation stopcock. L:270mm

Ref. D300 110 002



## Inner graduated operating sheath. 4.9mm

7Fr. operative canal. Rotatable irrigation stopcock. L:230mm

Ref. D300 110 001



## Outer graduated operating sheath. 5.5mm

Rotatable suction stopcock. Quick lock. L:216mm

Ref. D300 110 000



## Graduated pointed scissors. 7Fr.

Semi-rigid. Rotative. With flushing port  
L:370mm. Single action jaw

Ref. D300 120 065



## Graduated blunted scissors. 7Fr.

Semi-rigid. Rotative. With flushing port  
L:370mm. Single action jaw


Ref. D300 120 066



## Graduated biopsy forceps. 7Fr.

Semi-rigid. Rotative. With flushing port  
L:370mm. Double action jaw

Ref. D300 120 067

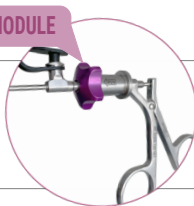

## Graduated forceps (alligator). 7Fr.

Semi-rigid. Rotative. With flushing port  
L:370mm. Double action jaw

Ref. D300 120 068



## Graduated IUD grasping forceps. 7Fr.

Semi-rigid. Rotative. With flushing port  
L:370mm. Double action jaw

Ref. D300 120 069



## Light cable

Autoclavable. Ø:3.5mm. L: 2.3m. Delivered with Storz connectors

Ref. D200 150 500



## EasyCare sterilization basket

Autoclavable. 540x240x70mm. With personalized silicone holders

Ref. D200 160 522



## Optional accessories:

### Bipolar electrode. Flexible. Button-type. 5 Fr.

Ref. D300 140 010



### Bipolar electrode. Flexible. Needle-type. 5 Fr.

Ref. D300 140 011

